

	Environmental Analysis Teaching and Research Laboratory	Date: 6/12/2022	SOP No. 23 v3
	Standard Operating Procedure	Title: IDEXX Quanti-Tray 2000 for E. coli and Total Coliform Bacteria Enumeration	
	Approved By: Marc Los Huertos	Revised: February 13, 2024	

1. Scope and Application

1.1 This standard operating procedure describes the test method for the collection and analysis of water samples for the enumeration of Escherichia coli (E. coli) and Total coliform bacteria.

1.2 The applications of this SOP are for...

2. Summary of Method

2.1 Surface water samples are collected in 120ml shrink-banded, sterile IDEXX bottles. An undiluted water sample will be analyzed from the sample collected. The Colilert reagent is added directly to the 100 ml undiluted sample. Both are mixed thoroughly to dissolve the reagent. The sample is transferred to QuantiTrays/2000 and sealed using the Quanti-Tray sealer. Samples are incubated at 35.0 ± 0.5 C for 24 hours. Results are reported as MPN/100mL.

2.2 The detection limit for this analysis is 1 Most Probable Number (MPN) per 100mL of sample.

Contents

1 Scope and Application	1
2 Summary of Method	1
3 Acknowledgements	3
4 Definitions	3
5 Biases and Interferences	3
6 Health and Safety	3
Safety and Personnel Protective Equipment	4
7 Personnel & Training Responsibilities	4

8 Required Materials and Apparati	4
9 Reagents and Standards	4
10 Estimated Time	4
11 Sample Collection, Preservation, and Storage	5
QA/QC	6
12 Procedure	6
Disposal	10
13 Data Analysis and Calculations	10
14 QC/QA Criteria	10
Data assessment and acceptable criteria	10
Corrective actions for out-of-control or unacceptable data	11
Disposal	11
15 Trouble Shooting	11
16 References	11

3. Acknowledgements

4. Definitions

4.1 Analytical batch: The set of samples processed at the same time

4.2 Control cultures: For each lot of medium, check analytical procedures by testing with known positive and negative control cultures. For example, E.coli is a positive control for this analysis and Staphylococcus aureus is a negative control.

4.3 Field duplicate (FD): Two samples taken at the same time and place under identical circumstances and that are treated identically throughout field and laboratory procedures. Analysis of field duplicates indicates the precision associated with sample collection, preservation, and storage as well as laboratory procedures.

4.4 Laboratory reagent blank (LRB): An aliquot of sterilized water treated as a sample in all aspects, except that it is not taken to the sampling site. The purpose is to determine if the analytes or interferences are present in the laboratory environment, the reagents, or the apparatus.

4.5 Laboratory duplicate (LD): Two aliquots of the same environmental sample treated identically throughout a laboratory analytical procedure. Analysis of laboratory duplicates indicates precision associated with laboratory procedures but not with sample collection, preservation or storage procedures.

5. Biases and Interferences

5.1 Water samples containing humic or other material may be colored. If there is background color, compare inoculated trays to a control tray containing only water (SM, 9223 A.)

6. Health and Safety

6.1 The analysis involves handling of freshwater samples that may contain live microorganisms and therefore pose some threat of infection. Laboratory personnel who are routinely exposed to such water samples are encouraged to protect themselves from water borne illnesses by wearing clean disposable gloves and washing their hands frequently.

6.2 The Colilert reagent is not hazardous according to the manufacturers material safety data sheet. The manufacturer does recommend wearing gloves and safety glasses while using this reagent and washing hands after use.

Safety and Personnel Protective Equipment

Water quality can range from high quality and even potable to very low with assorted pathogens or even toxicity. Thus, it's important to take every precaution possible to protect yourself.

7. Personnel & Training Responsibilities

7.1 Laboratory and field personnel shall have a working knowledge of this analytical procedure and will have received training from a faculty, staff, or student knowledgeable of the proper sample analysis procedures.

7.2 Researchers training is required before this the procedures in this method can be used.

7.3 Researchers using this SOP should be trained for the following SOPs:

- SOP 01 Laboratory Safety
- SOP 02 Field Safety

8. Required Materials and Apparati

8.1 10.1 Sterile, shrink-wrapped 100ml IDEXX bottles- Catalog NO: WV120SBT-200

8.2 Quanti-Tray Sealer:IDEXX Laboratories, Inc., Westbrook, ME -Serial NO: QT1757, Model NO: B9-10894-04 .

8.3 Incubator- Catalog NO: 150152632, Serial NO: IMC40490982

9. Reagents and Standards

9.1 Colilert reagent- Catalog NO: WP100l-18

9.2 ?? Quanti-Tray/2000: 100 trays containing 97 wells each IDEXX Laboratories, Inc., Westbrook, ME -Catalog NO: WQT2K

10. Estimated Time

10.1 This procedure requires XX minutes...

11. Sample Collection, Preservation, and Storage

- 11.1** Arrive at site and record site number, date and time.
- 11.2** Immerse the thermometer or YSI handheld in the water and leave immersed five minutes before reading temperature. Avoid disturbing the bottom with the thermometer at the sample site.
- 11.3** Label bottle with location (geographic area name and stream or lake name), date, time, water temperature and sampler's initials. Label bottle before immersion using a black permanent marker or pre-printed labels. QVIR Bacteria Lab, State Certified Lab, purchases only certified sterile, 100 ml, sealed containers from IDEXX.
- 11.4** Use latex gloves when handling bottles during sampling. Fingers contain contaminants such as nitrates. Bug repellents or sunscreen are particularly troublesome as contaminants. Once the gloves are on, be careful not to touch your face, the ground, or anything but the bottles.
- 11.5** The sample should be taken from flowing, not stagnant water, facing upstream.
- 11.6** Be sure to immerse the bottle completely, 10 cm (4 inches) deep, with mouth of bottle pointing upstream, so no water flows over your hand into the bottle. Be sure the bottle does not get near the bottom of the stream where sediments can be disturbed. Water samples should be collected 6-12 inches below the water surface. Fill bottle, to the 100ml line indicated, on first immersion, pour off the excess and cap. Do not under fill or over fill, do not redunk. If too much water is poured off, redo sample with new 100 ml container.
- 11.7** Do not touch bottle mouth or inside of cap. Be careful not to contaminate the sample with surface film, contact with human skin, breathing in/on the bottle or cap, etc. If stream is too shallow to immerse bottle fully, collect as much as possible, being very careful not to touch the bottom. Note depth on field notes.
- 11.8** Collect one "duplicate" sample every two weeks (sampling frequency). Sample sites chosen for duplicate sampling are selected at random among sites sampled. When a duplicate sample is selected for the site, repeat procedures as with normal stream samples. The duplicate is the second sample when two samples are collected. Duplicates document repeatability of individual sample collections and reproducibility of laboratory results.
- 11.9** Samples are analyzed in the Biogeochemistry Lab. Keep samples cool while transporting. Store at 4 C until analysis, but do not freeze. The maximum holding time is 6 hours.
- 11.10** For each sample, the location number, bottle numbers used and time collected will be recorded in the field sample log.

11.11 The samples will be kept in the possession of researcher personnel who both collect and analyze the samples.

QA/QC

11.12 Accuracy: Initial analyst demonstration of capability and for each new lot of Quanti-Tray/2000, analyze the following:

11.13 Check each new lot of Colilert. Shine the ultraviolet lamp on the media snap packs. If the lot is fluorescent it will be discarded.

11.14 Dissolve one packet in 100 ml distilled water. Do not incubate. Check for fluorescents.

11.15 Analyze sterile reagent water blank with each batch of samples to verify that there is a negative result from 24-28 hours

11.16 Gravimetrically check each new lot of sterile, transparent, nonfluorescing 100-ml vessels to ensure the 100-mL fill line is accurately represented on the vessel.

12. Procedure

Good laboratory practices will be used for this procedure.

12.1 Turn on the Quanti-tray sealer.

12.2 Set the incubator to 35 °C.

12.3 Add Colilert reagent to 100 ml water sample.

12.4 Place the lid back on the sample and shake well till the reagent has dissolved.

12.5 Use a permanent felt tip marker to label your Quanti-tray with information identifying your sample.

12.6 Open the Quanti-tray by holding it at the top, with the well side facing your palm. Push the top edge of the Quanti-tray in with your free hand while squeezing the tray into a circle. Open the tray by gently pulling the foil tab away from the plastic side. Be careful not to tear the tab. Do not touch the inside of the tray.



Figure 1: Adding Reagent

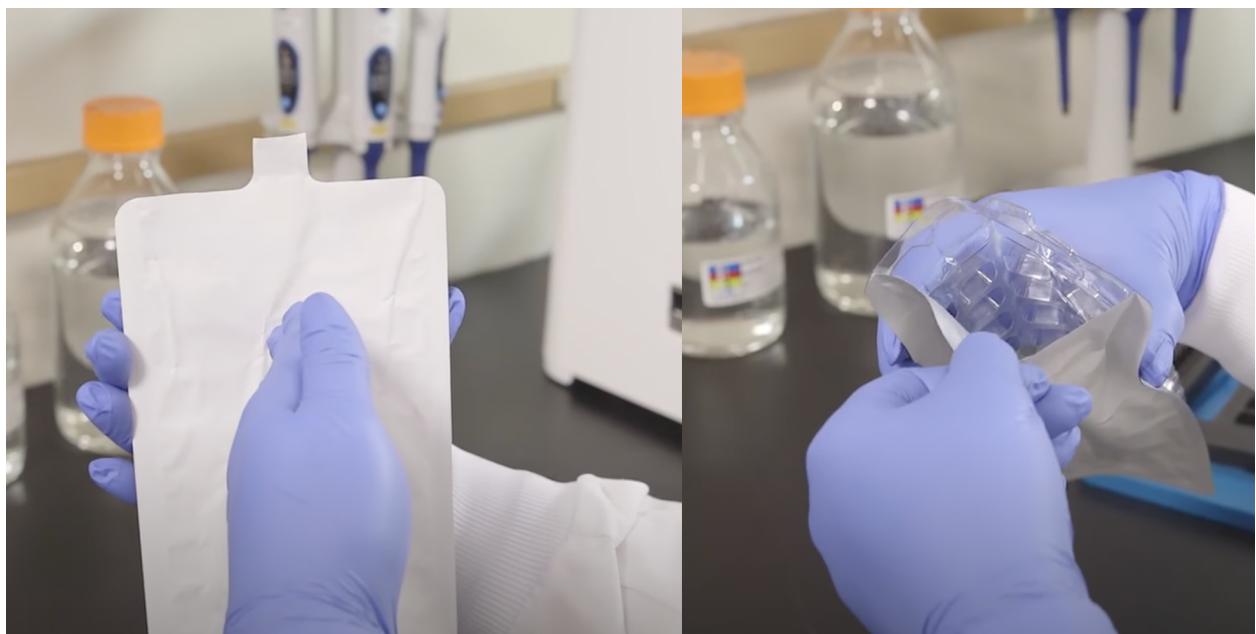


Figure 2: Opening and Pouring Sample into Quanti-Tray



Figure 3: Pouring Sampling into Quanti-Tray

- 12.7** Pour the sample mixture into the Quanti-tray, avoiding contact with the foil tab.
- 12.8** Place the filled Quanti-tray on the rubber insert, making sure each well fits its corresponding hole.
- 12.9** Push the insert with the tray into the sealer, non foil tab end first, until the sealer grabs the tray and pulls it into the slot. If you need to reverse the motor, press and hold the reverse button. However, do not reverse the motor if the rubber insert is completely inside the sealer. The sealer distributes the sample mixture into the Quanti-tray wells, seals the wells, and then partially ejects the sealed tray.



Figure 4: Pressing the Tray into the Rubber Insert

- 12.10** Remove the rubber insert and tray from the sealer.

12.11 Incubate the sealed tray for 24 hours at 35 degrees C plus or minus 0.5 degrees C.

12.12 After incubation, mark each yellow well with a permanent marker, including the large overflow well at the top. To determine the number of total coliforms, count the wells that are equal to or greater than the comparator, and then refer to the MPN.



Figure 5: Counting Coliform and E Coli Wells (Be careful to avoid eye damage with the use of the UV lamp.)

12.13 To determine the number of E Coli, view the Quanta-tray under a UV lamp. Count the fluorescent wells that are equal to, or greater than the comparator, and then refer to the MPN table.

Table 1: need to add caption

Organism	Result
Escherichia coli	Yellow wells, fluorescence
Klebsiella pneumoniae	Yellow wells, no fluorescence
Pseudomonas aeruginosa	Clear wells, no fluorescence
Method Blank	Clear wells, no fluorescence

12.14 A count of 110 MPN/100 ml is regarded as low risk; 11100 MPN/100 ml is medium risk. Finally, an E. coli count greater than 100 MPN/100 ml is judged high risk.

12.15 Report results on the bench sheet.

12.16 The completed bench sheet should be reviewed by the researcher, the laboratory director and the QA manager.

Disposal

12.17 All materials must be autoclaved prior to disposal and workspaces thoroughly disinfected.

12.18 Dispose of media in accordance with Good Laboratory Practices.

13. Data Analysis and Calculations

13.1 For each sample analyzed, including quality control samples, record the number of small and large positive wells and the MPN in the appropriate places on the bench sheet (see below). Calculate precision for duplicate analyses using equation 1.

13.2 Equation 1.

$$Precision(asRPD) = \frac{(AB) * 100}{(A + B)/2} \quad (1)$$

Where: A = MPN from aliquot A and B = MPN from aliquot B

13.3 The following results should be observed:

14. QC/QA Criteria

Data assessment and acceptable criteria

14.1 The analyst should review all data for correctness (e.g., use of MPN table).

14.2 Record the precision values as RPD on the bench sheet.

14.3 The desired precision is 20

14.4 The desired detection limit is 1 MPN/100mL

14.5 The completed bench sheet is reviewed by the analyst's supervisor or the QVIR Lab Director

Corrective actions for out-of-control or unacceptable data

14.6 The results for precision and blank data are compared to the acceptable values for this analysis; 20% and 1 MPN/100mL, respectively.

14.7 If a precision value exceeds 20% then the analyst should write in the comments section of the bench sheet: These data are associated with an out-of-control duplicate analysis. The UCL = 20%. Note: UCL is the Upper Control Limit (i.e., 20

14.8 If a blank value exceeds 1 MPN/100mL then the analyst should write in the comments section of the bench sheet: These data are associated with a blank value that exceeds the detection limit of 1 MPN/100mL.

14.9 The samples cannot be reanalyzed because the sample volume will be depleted after the initial analysis.

14.10 If data are unacceptable for any reason, the analyst should review their analytical technique prior to conducting this analysis again.

Disposal

14.11 All materials must be autoclaved prior to disposal and workspaces thoroughly disinfected.

14.12 Dispose of media in accordance with Good Laboratory Practices.

15. Trouble Shooting

TBD!!!

16. References

16.1 APHA, AWWA, WEF. (2012) Standard Methods for examination of water and wastewater. 22nd American Public Health Association (Eds.). Washington. 1360 pp. (2014).